


CT ACR 464 Phantom Body Ring & Extensions

Expand the functionality of your CT ACR 464 Phantom.

- Simulate larger patients
- Represent scatter effects from widebeam CT scanners

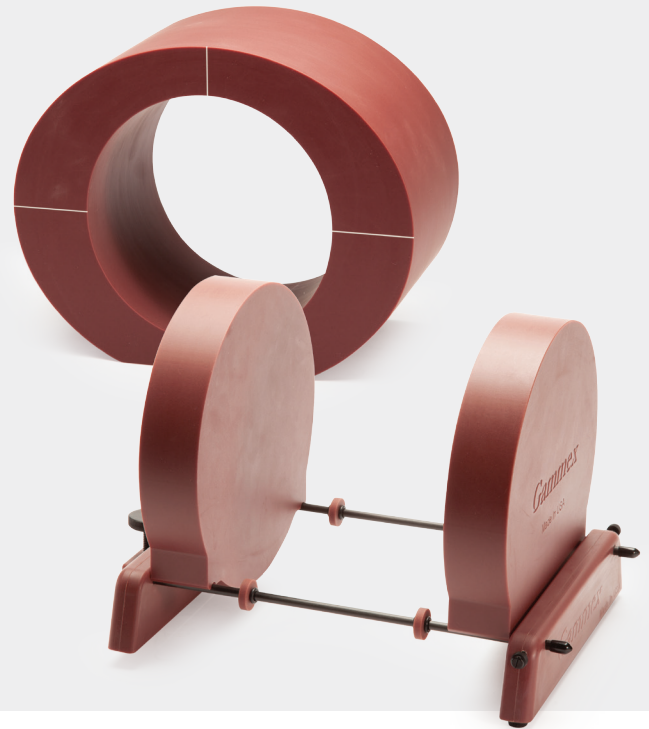

CT ACR 464 Phantom Body Ring

The CT ACR 464 Phantom can easily be inserted into the Body Ring adapter to permit fast imaging QA for a better indication assesment of a scanner's performance with large patients.

Benefits

- Expand the use of the CT ACR 464 Phantom far beyond accreditation
- Fast setup and easy to use
- Mimics a human torso to better represent the typical scanner imaging area and beam hardening effects

Specifications

Material:	CT Solid Water®
Diameter:	33.0 cm (13 in)
Length:	26.4 cm (10.4 in)
Depth:	16.0 cm (6.3 in)
Weight:	6.1 kg (13.5 lbs)

CT ACR 464 Phantom Extensions

The Extension plates accurately represent scatter effects from widebeam CT scanners thus eliminating the need to add a water bolus or other material when measuring the CT ACR 464 .

Benefits

- Kit includes 2 extension plates and stand
- Fast setup and easy to use
- Custom cases available for safe storage and easy transport

Specifications

Material:	CT Solid Water®
Extension Plate Diameter:	20.0 cm (7.9 in)
Extension Plate Thickness:	4.0 cm (1.6 in)
Overall length of 464 Phantom plus Extension Plates:	24 cm (9.4 in)
Weight:	3.2 kg (7.0 lbs)

